



**Para Mobility**

Making everyday living easier



 horizon.

# Montar Wall Mounted Change Table

Product Code: 7-MONTAR

---

## User Manual and Monthly Maintenance Checklist

**IMPORTANT NOTICE:** This document must be forwarded to your architect, builder or installer to ensure the site is prepared correctly, as any misuse or understanding of the installation of this product could be costly, cause accidents or be fatal.

Please carefully read instructions before using this change table.

## CONTENTS

User Manual  
Monthly Maintenance Manual  
Monthly Checklist  
Warranty

## INTRODUCTION

This instruction manual contains information on the instructions for use and explanations of the signs and symbols used. These instructions for use contain instructions for operating the medical care bed.

These instructions for use may contain inaccuracies or printing errors. The information provided here will be updated periodically and changes due to product maintenance will be included in later editions. Changes or improvements are possible at any time without prior notice. If you have any questions, please contact our customer service team.

The instructions for use must be read and used by every person who operates the medical bed.

## PURPOSE

The intended purpose of the Horizon Montar Wall Mounted Change Table is an electric Examination couch/table for examination and/or treatment that may be fitted with some basic functions, e.g. raise & lower, and is used in, e.g. examination rooms or doctors' surgeries. These are electrically powered to allow the bed to be placed in various height positions. It is a medical device intended for human medical use only.

The Horizon Montar Wall Mounted Change Table features adjustable motorised height settings.

The intended user of the Horizon Montar Wall Mounted Change Table is any individual receiving regular changing or being examined.

The bed should only be used by individuals trained or experienced in proper handling. Users must be instructed on its correct use and familiarize themselves with these instructions. Common users include nursing staff, physicians, therapists, family members, technicians, and cleaning staff.

## Symbols Explained:


---



Manufacturer



## PRODUCT SERIAL NUMBER PLATE

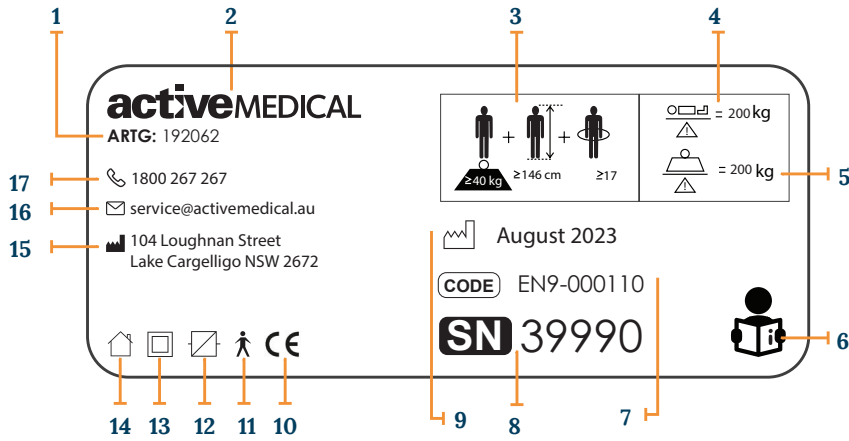

Serial Number Plate

### Serial Number Plate Explained:

|    |   |
|----|---|
| 1  | ARTG Number   |
| 2  | Company Logo  |
| 3  | Adult Identification (minimum weight, minimum height, minimum BMI)                                    |
| 4  | Maximum Occupant Weight   |
| 5  | Maximum Safe Working Load   |
| 6  | Read Instruction Manual - Follow Instructions For Use   |
| 7  | Product Code  |
| 8  | Product Serial Number   |
| 9  | Manufacture Date  |
| 10 | CE Marking  |
| 11 | Type B Application Part (classification according to the degree of protection against electric shock) |
| 12 | Switching Power Supply Unit   |
| 13 | Device Protection Class II (classification according to the type protection against electric shock)   |
| 14 | For Indoor Use Only   |
| 15 | Manufacturer Address  |
| 16 | Company Service Email Address   |
| 17 | Company Phone Number  |



## IMPORTANT

**WARNING:** The word warning refers to a hazard to unsafe practice that may cause severe injury or death to you or another person.


## NOTES

1. Horizon shall not be liable for loss, cost, damages or expenses incurred by the client or any other person or company, resulting from the use of any information or interpretation given in this report. In no case shall Horizon be liable for consequential damages, including but not limited to, lost profit, damages for failure to meet deadlines and lost production arising from this report. This document shall not be reproduced except in full and relates only to this product.
2. It remains the responsibility of the client to ensure that this product tested is representative of the entire product batch.
3. Horizon shall take no responsibility for the procurement and authenticity of the test product as described herein.
4. This report is specific to the test items in their state at the time of testing. It should not be taken as a statement that all products in all states of repair would also perform in the same manner.
5. Horizon shall take no responsibility for the interpretation of the procedures or calculation methods as provided herein or for the appropriateness or validity of the test procedures for the test items described and reported.
6. Horizon shall take no responsibility for the installation procedures used for the test items as described herein.
7. Horizon shall take no responsibility for the Safe Working Load limits of 2000N and 4000N as specified in AS/NZS 3200 Part 2.38-2007 (Electrically and Manually Operated Medical Beds for Adult Use), Section 21.3.101 Mechanical Strength & Section 28.4 Suspended Masses respectively.


**WARNING:** No patient, service personnel or any person must lie underneath the Montar while power is turned on. If not obeyed, this may result in serious bodily injury or death if the bed controls are actuated.

**WARNING:** No person should ever be left on the change table unattended.

## Montar Specifications:

|  |   |               |
|--|---|---------------|
| Overall Width Open/Closed  | 350/1030 mm   | 350/1030 mm   |
| Overall Length   | 1400 mm   | 1880 mm       |
| Sleeping Surface   | 700 x 1400 mm   | 700 x 1880 mm |
| Product Weight   | 100 kg  | 100 kg        |
| SWL  | 200 kg  | 200 kg        |
| Low/High Position <b>**approx**</b>  | 385/985 mm  | 385/985 mm    |
| TGA Number   | ARTG 135743   |               |
| Manufacturer  | Active Medical, 104 Loughnan Street, Lake Cargelligo NSW 2672 |               |

## MONTAR POSITIONS

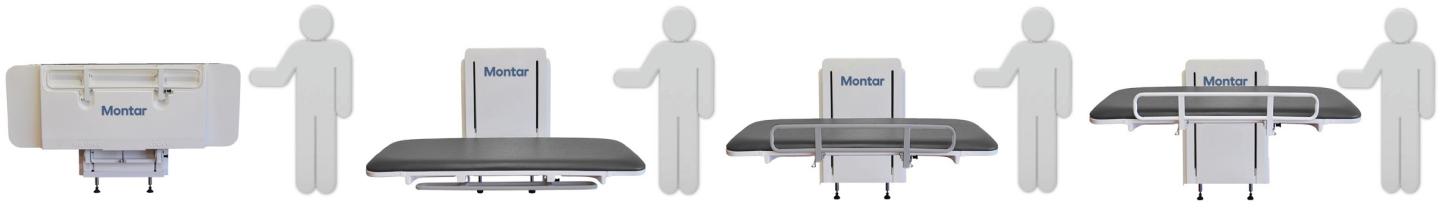

## EXPLANATION OF FEATURES



### **Space Saving Position**

The Montar in its space friendly position. When it's not in use, it can easily fold to save space and avoid trip hazards.



### **Unfolding Bed Base**

**LOCK** - to fold the Montar's bed base to its lying position, put one hand on the top of the bed base - unlock - and slowly lower.



### **Handset Functions**

Handset showing up/down buttons for raising or lowering the padded surface.



### **Patient Entering Position**

The Montar in it's patient entering position - this is the lowest position with the rail tucked underneath for easy access.

**WARNING:** No person should be left on the change table unattended.





### **Moving Rails to Fixed Position**

Engaging rail to its upward fixed position by simply pulling the rail towards you as well as upwards. This will lock automatically once the rail is 90 degrees to the bed base.



### **Lowest Position**

The Montar positioned in its lowest position with rail in its lowest position.

**WARNING:** Ensure safety rail is secured whenever a patient is lying on the surface.



### **Unlocking Rail**

To unlock the rail, pull the lock knob, gently lower and push the rail back underneath the bed base to its home position.



### **Highest Position**

The Montar in its highest position.

**WARNING:** DO NOT enter or leave the Montar whilst it's in the highest position.



## MONTHLY MAINTENANCE MANUAL AND CHECKLIST (PARTS TO CHECK ON THE MONTAR)



### **Wear and Tear**

Look for wear and tear on the walls, cracking cornices and listen for creaking (with weight on the Montar).



### **Structural Check**

Ensure the Montar is fastened securely to the wall and fasteners are tight. Ensure there are no cracks in the walls or framework of the Montar.



### **Safety Rail Check**

Ensure the safety rail bolts are secure and the rail is functioning properly.



### **Bed Base Check**

Check the bed base hinge bolts are secure and tight.



### **Gas Strut Check**


Ensure gas strut is secure, tight and functioning properly.



### **Feet Check**

Make sure the feet are firm and there is weight onto the ground - if the feet are loose, they will need rescrewing downwards to force weight back to the ground.



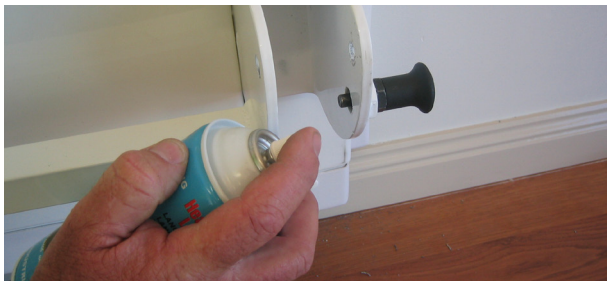

### Moving Parts Check

Spray moving parts and locks with some lubricant spray - this will enhance longer life to wearing parts.



### Oil Guide Check

Spray and oil guide tubes through the cutouts on the wall unit cover.

## CLEANING THE MONTAR WALL MOUNTED CHANGE TABLE

- Cleaning can be done often - using domestic detergents.
- Do not use acids, alkaline or any other chemical solvents.
- When cleaning, be careful of entrapment, as it could cause physical injury.

## MONTHLY CHECKLIST

Please use the below checklist when completing your monthly checks.

### Record of Performed Service Procedures and Regular Service Maintenance

|               | JAN | FEB | MAR | APR | MAY | JUN | JULY | AUG | SEPT | OCT | NOV | DEC |
|---------------|-----|-----|-----|-----|-----|-----|------|-----|------|-----|-----|-----|
| Wear & Tear   |     |     |     |     |     |     |      |     |      |     |     |     |
| Structural    |     |     |     |     |     |     |      |     |      |     |     |     |
| Safety Rail   |     |     |     |     |     |     |      |     |      |     |     |     |
| Bed Base      |     |     |     |     |     |     |      |     |      |     |     |     |
| Gas Strut     |     |     |     |     |     |     |      |     |      |     |     |     |
| Feet          |     |     |     |     |     |     |      |     |      |     |     |     |
| Moving Parts  |     |     |     |     |     |     |      |     |      |     |     |     |
| Oil Guide     |     |     |     |     |     |     |      |     |      |     |     |     |
| Battery Check |     |     |     |     |     |     |      |     |      |     |     |     |

**HAVE YOU TICKED EVERY MONTH THIS YEAR?**







## HORIZON MONTAR WALL MOUNTED CHANGE TABLE WARRANTY

### Conditions

From the date of purchase, your Horizon Wall Mount Change Table is guaranteed for 10 years on the frame, 5 years on the electrics and 1 year on the battery.

Professional service of product to be carried out annually and documented in the Record of Performed Service Procedures and Regular Service Maintenance.

This warranty applies only on the parts manufactured or approved by Active Medical.

### Exclusions

This warranty does not extent those items that may need replacement due to normal wear and tear, or to damaged caused by misuse or accident, for which Active Medical or it's distributors cannot be held responsible.

If your change table should need attention due to failure as defined under the warranty above, please contact the distributor from whom you purchased your bed.

Please fill in the below and retain for your own records.

|                                   |
|-----------------------------------|
| <b>Date Purchased</b>             |
| <b>Model</b>                      |
| <b>Serial Number(s)</b>           |
| <b>Distributor</b>                |
| <b>Distributor Address</b>        |
| <b>Distributor Contact Number</b> |





✉ [hello@horizoncare.au](mailto:hello@horizoncare.au)

☎ 1800 267 267