


# PLATYPUS POOL STEPS

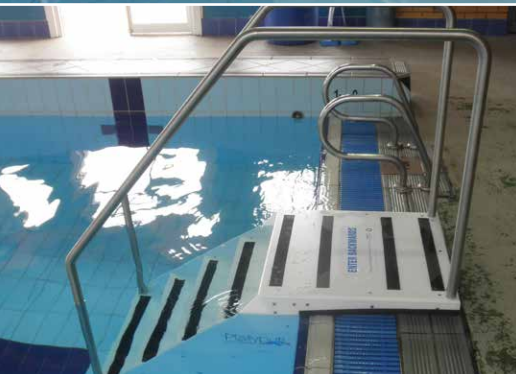

## Open your pool to everyone

Swimmers of all types benefit from the innovative design of the Platypus Pool Steps - elderly, arthritic, expectant mothers, young learners or anyone who is uncomfortable with, or has difficulty using traditional vertical swimming pool ladders.

The Platypus Pool Steps offer no deck modifications or alterations to install. The Platypus Pool Steps can be placed in and out of your pool as many times as needed - safely and easily by its removable dolly wheels.

Manufactured in Sydney, using stainless steel hand rails and corrosion resistant HDPE plastic treads. The Platypus Pool Steps are your long term solution to your swimming pool access problems.



## PLATYPUS POOL STEPS **PRODUCT SNAPSHOT**

### FEATURES

- Makes your pool more accessible
- Custom made to fit your pool perfectly
- Can be fitted to almost any in-ground pool or decked above ground pools
- Easy to pull out of the pool and place back in when needed
- Can be installed in the corner or on the side of your pool
- SWL 150kg



Phone (02) 9651 4446

Delivery, installation + servicing Australia wide  
[sales@paramobility.com.au](mailto:sales@paramobility.com.au)

**[www.paramobility.com.au](http://www.paramobility.com.au)**



Disability Discrimination Act 1992

Proud supporters of







**Earthing:** As per AS3000 of the wiring rules standards – “In certain circumstances” it may be necessary to have all conductive metal objects equally potentially bonded (earthed), this includes Pool Hoists.


**We recommend you seek the advice and opinion of a qualified electrician.**

## TOOLS REQUIRED


- Two spanners

## INSTALLATION INSTRUCTIONS PELICAN POOL STEPS

- Your Platypus Pool Steps will arrive on a pallet.
- Carefully unpack the steps and transfer using the included dolly wheels to the area set out for the steps.
- Lower the steps into the pool, this will require the help of someone in the pool.
- All four adjustable feet must be touching the bottom of the pool.
- They can be adjusted by using two spanners.
- The top platform of the steps that sits on the pool deck must also be sitting level on the pool deck.

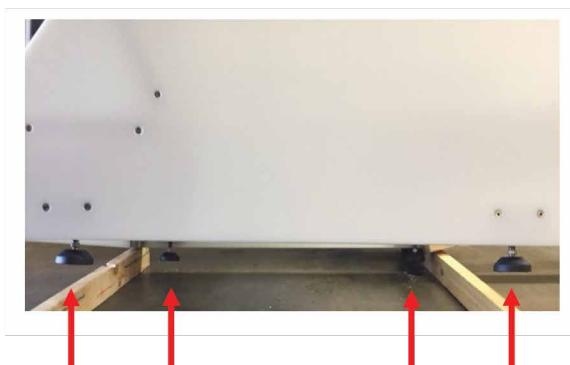

All four feet can be adjusted by using two spanners.

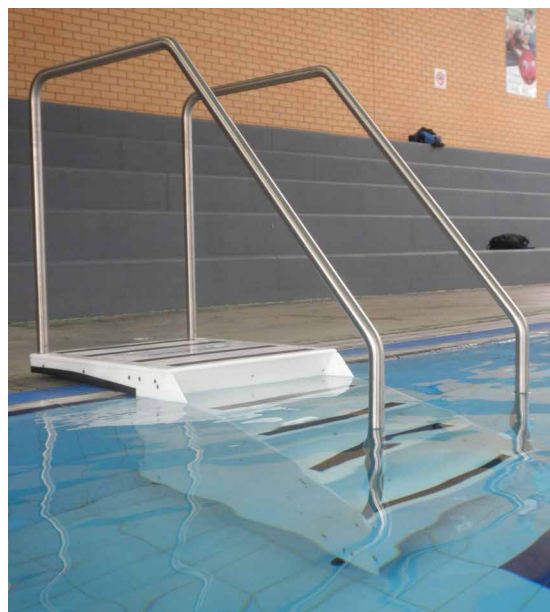



## USER INSTRUCTIONS PELICAN POOL STEPS

### ENTERING THE POOL

The Platypus Pool Steps have been designed to be used as a ladder. It is critical that all users are instructed to enter and exit the pool facing the same direction – away from the pool.

Never allow the user to enter or exit the pool facing forwards. Serious injury may occur.

### CORRECT

fig. 3



fig. 4

### INCORRECT



### REMOVING THE STEPS

Removing the steps is always a two-person job.

Attach the Dolly Wheels to the standing platform before moving. Please use caution and care when removing the steps from the pool so as not to injure yourself.







## CLEANING & MAINTENANCE

Para Mobility only use 316 grade Stainless Steel to manufacture our products, due to the corrosion or oxidation (rust) resistant properties.

Tea Staining is a discolouration of the surface of the stainless steel, and occurs in higher temperatures and humid conditions. Whilst having the same look as rust, Tea Staining does not affect the structural integrity or the longevity of the material. Visually it is a discolouration of the surface, which tends to follow the "grain" of the surface finish.

To ensure your pool steps continue to look great and avoid tea staining we recommend washing the area regularly using clean fresh water and mild detergent using a cloth or soft brush.

After washing rinse off with clean water and wipe the surface dry with a soft absorbent cloth.

**DO NOT** use abrasive cleaners or a scouring pad of any type as they may cause scratching, damage, excessive wear and dull the finish.


**DO NOT** use waxes or furniture polish.

## USER INSTRUCTIONS PELICAN POOL STEPS

### TRANSPORTING THE STEPS

Never drag the steps across the pool deck.


Attach the Dolly Wheels to the standing platform and hold underneath the bottom step when moving the steps.

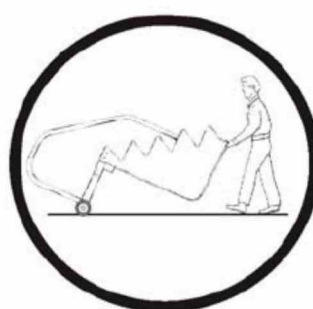

**CORRECT**



**NOT CORRECT**

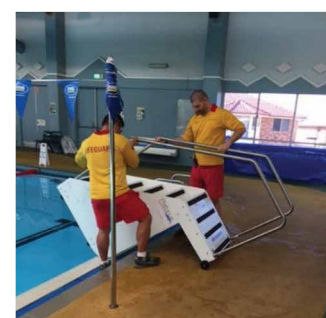

## User Instructions

### Entering the Pool

The Platypus Pool Steps have been designed to be used as a ladder. It is critical that all users are instructed to enter and exit the pool facing the same direction – away from the pool.

Never allow the user to enter or exit the pool facing forwards. Serious injury may occur.


**CORRECT**

fig. 3



fig. 4

**INCORRECT**

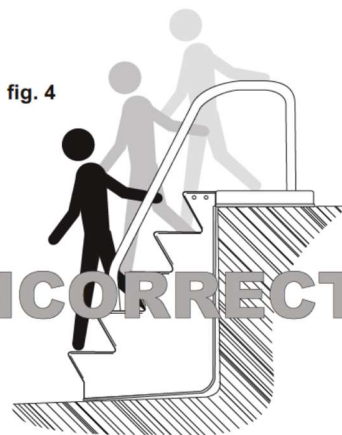

### Removing the Steps

Removing the steps is always a two-person job.


Attach the Dolly Wheels to the standing platform before moving.

Please use caution and care when removing the steps from the pool so as not to injure yourself.



### Transporting the Steps

Never drag the steps across the pool deck. Attach the Dolly Wheels to the standing platform and hold underneath the bottom step when moving the steps.

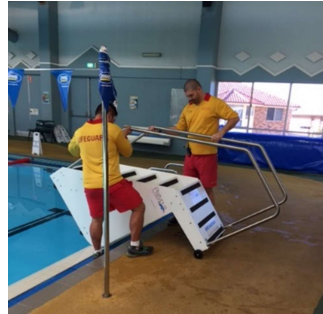

**CORRECT**



**NOT CORRECT**

### Cleaning and Maintenance

Only use a non-abrasive cleaner when cleaning the Platypus Pool Steps to avoid scratching and damage to the finish. Never use waxes or furniture polish.

**Manufacturers and Importers of Specialised Disability Equipment since 1985**

**A:** 31 / 276 New Line Road, Dural NSW 2158 **P:** 02 96514446 **E:** sales@paramobility.com.au

**W:** www.paramobility.com.au