

Pelican Pool & Spa Hoist



Designed and manufactured in Australia by Para Mobility, and tested in accordance with ASA 3581, the Pelican Pool & Spa Hoist allows a person with any level of disability to have access to both above and below ground spas, swimming & hydrotherapy pools – without any manual lifting. It is suitable for use in public pools, leisure centres, private homes, nursing homes and schools.

The unit can be mobile or fixed permanently at the edge of the pool or spa and is easily installed on the pool deck. Our team work with you to ensure the Pelican Pool and Spa Hoist is made to meet your pool needs.



FEATURES

- TGA approved
- NDIS confirmed
- Easily wheeled from one area to another for multiple use and or storage (mobile unit)
- No manual lifting required
- Press button hand switch operation
- Standard lifting capacity 150kgs and 200kg
- Cover included
- Easily installed
- Complies with Disability Discrimination Act (Access to Premises – Buildings) under the Building Code of Australia

SEATING OPTIONS


- Floating sling made by Para Mobility
- Solid plastic seat for clients with upper body control
- Pelican Aquatic Wheelchair – designed and manufactured by Para Mobility, a great option as it can be used around the pool area
- Platform lift (Meets Access to Premises Standards 2010 for pools greater than 70m perimeter)

OPTIONAL EXTRAS

- Power rotation – push button operation for independent use
- Earthing ability
- Bluetooth handset

SPECIFICATIONS

- Safe Working Load: 150kg and 200kg
- Batteries: Linak – 24v re-chargeable battery
- Rotation: Manual – 360 degrees. Power rotation available.
- Total Weight: From 60kg
- Floor Mounting Plate: 450 diameter
- Construction: Stainless Steel, powder coated to colour of choice
- Hand control: Linak – push button
- Installation: Anchored to pool deck with Stainless Steel anchors and eye bolts

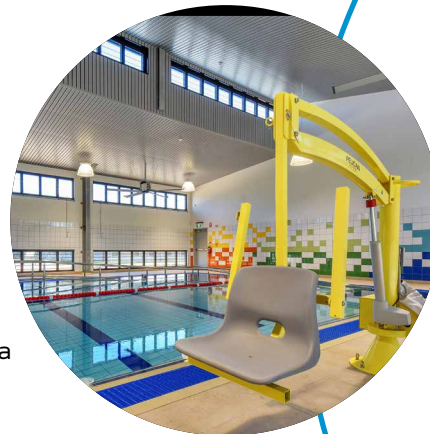

Freedom to explore

